

Auto report file naming

QC

- Use CV_export_app to export the table
- Filename structure:

QC_PTK.csv



File type. Should be QC.







Assay type. Options: PTK or STK

Phosphosite Analysis (Limma)

- Running a stats app in Tercen will automatically export the table needed for reporting. Rename this file.
- PTK and STK files have to contain the exact same supergroup name!
- Filename structure:

Limma_PTK_01_Supergroup.csv

-  Stats type: Limma
-  Assay type. Options: PTK or STK
-  Order, defines ordering in report.
-  Supergroup: the name of the supergroup **factor**.
NOT the name of the supergroup! This is NOT displayed in Table 1.
Must be unique. E.g. Supergroup or Test condition.
Both supergroup names and comparisons come from the file.





Output:

Assay Type	PTK	
Comparisons	Up	Down
1-A219 - T1 vs Control ^a	8	11
1-A219 - T2 vs Control ^a	0	7
1-A219 - T2 vs T1 ^a	1	1
^a Significance was obtained using Limma, $p < 0.05$.		

Phosphosite Analysis (T-Test supergroup)

- Running a stats app in Tercen will automatically export the table needed for reporting. Rename this file.
- PTK and STK files have to contain the exact same supergroup name!
- Filename structure:

TT_PTK_01_T vs C.csv

-  Stats type: TT
-  Assay type. Options: PTK or STK
-  Order, defines ordering in report.
-  Comparison, e.g. T vs C.
Comparison is displayed in Table 1.
Supergroups come from the file.

Output:




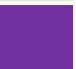
Table 1: Phosphosite Analysis

Assay Type	PTK	
	Up	Down
Comparisons		
Sgroup1 - T vs C ^a	23	0
Sgroup2 - T vs C ^a	4	5
^a Significance was obtained using a two-sided unpaired Student's T-test, $p < 0.05$.		

Phosphosite Analysis MTvC

- MTvC filenames need to be different if you want to include multiple MTvC files in one report
- Running a stats app in Tercen will automatically export the table needed for reporting. Rename this file.
- Filename structure:

MTvC_PTK_01_10h.csv

-  Stats type: MTvC
-  Assay type. Options: PTK or STK
-  Order, defines ordering in table.
-  Name of Supergroup. E.g., timepoint, cell line. This is displayed in Table 1. Comparisons come from the file.

Output:

Table 1: Phosphosite Analysis

Assay Type Comparisons	PTK	
	Up	Down
10h - T1 vs <u>Control</u> ^a	0	32
10h - T2 vs <u>Control</u> ^a	0	8
10h - T3 vs <u>Control</u> ^a	10	0

^aSignificance was obtained using a one-way ANOVA followed by a post-hoc Dunnett's test, $p < 0.05$

Kinase Analysis (UKA 2023)

UKA MTvC filename structure:

UKA_PTK_01_ukam-Sgroup1 vs C.csv

↓ Autoreport generates „old” data format with these names:

UKA_PTK_01_Sgroup1 - T1 vs C.csv

UKA_PTK_02_Sgroup1 - T2 vs C.csv

- File type: UKA
- Assay type (PTK or STK)
- Order. Defines ordering in table.
- Uka types (ukam or ukat)
- Supergroup. If you don't have supergroup, give it a neutral name, e.g. Test
- Comparison. In format XvsY or X vs Y.

UKA TGC filename structure:

UKA_PTK_01_ukat-T vs C.csv



UKA_PTK_01_Sgroup1 - T vs C.csv

UKA_PTK_02_Sgroup2 - T vs C.csv

UKA_PTK_01_ukat-Sgroup2 vs Sgroup1.csv



UKA_PTK_01_Control - Sgroup2 vs Sgroup1.csv





UKA_PTK_02_Test - Sgroup2 vs Sgroup1.csv

Comes from the file/Test Condition, not from the filename.

Kinase Analysis (UKA 2022)

- Export report summary from UKA
- Filename structure:





UKA_PTK_01_T1vsT2.csv

-  File type: UKA
-  Assay type. Options: PTK or STK
-  Order. Defines ordering in table.
-  Comparison. In format XvsY or X vs Y.

Phosphosite Analysis (single MTvC or TT) deprecated!

- Running a stats app in Tercen will automatically export the table needed for reporting. Rename this file.
- Filename structure:

MTvC_PTK.csv

-  Stats type. Options: MTvC or TT
-  Assay type. Options: PTK or STK
-  Order, only for TT. Defines ordering in table.
-  Comparison, only for TT. In format XvsY or X vs Y.

QC

- Export **flat file** from BioNavigator after log2 transformation
 - Array factors: Barcode, Array, and Test condition
 - Spot factor: ID
 - Quantitation type: S100-logTransformed
- Filename structure:

QC_PTK.txt



File type. Should be QC.








Assay type. Options: PTK or STK

Phosphosite Analysis (T-Test supergroup)

- Take tab-delimited text files from BioNavigator of LogFC and P values from the DataExports folder
- PTK and STK files have to contain the exact same supergroup name!
- Filename structure:






TT_PTK_01_Group_LogFC.txt

-  Stats type. Options: MTvC or TT
-  Assay type. Options: PTK or STK
-  Order, only for TT. Defines ordering in report.
-  Group, defines the name of the supergroup. Must be unique. Comparisons come from annotation.
-  File type. Options: LogFC.txt or p.txt

Phosphosite Analysis multiple MTvCs

- MTvC filenames need to be different if you want to include multiple MTvC files in one report
- Take tab-delimited text files from BioNavigator of LogFC and P values from the DataExports folder
- Filename structure:





MTvC_PTK_01_10h_LogFC.txt

-  Stats type. Options: MTvC
-  Assay type. Options: PTK or STK
-  Order, only for TT. Defines ordering in table.
-  Group name, indicates the type of MTvC. E.g., timepoint, cell line.
-  File type. Options: LogFC.txt or p.txt

Kinase Analysis

- Export report summary from UKA
- Filename structure:






UKA_PTK_01_T1vsT2.txt

-  File type: UKA
-  Assay type. Options: PTK or STK
-  Order. Defines ordering in table.
-  Comparison. In format XvsY or X vs Y.

Phosphosite Analysis (single MTvC or TT) deprecated!

- Take tab-delimited text files from BioNavigator of LogFC and P values from the DataExports folder
- Filename structure:

MTvC_PTK_LogFC.txt

-  Stats type. Options: MTvC or TT
-  Assay type. Options: PTK or STK
-  Order, only for TT. Defines ordering in table.
-  Comparison, only for TT. In format XvsY or X vs Y.
-  File type. Options: LogFC.txt or p.txt